




6197 Cornerstone Court, Ste 102
San Diego, CA 92121

Product Specifications

UltraPolymer Universal Anti-Mouse/Rabbit IgG (H&L) – HRP

Description:	UltraPolymer Goat anti-Mouse/Rabbit IgG (H&L) conjugated to HRP: Goat anti-Rabbit IgG (H&L) conjugated to HRP, affinity purified, min x w/ human, mouse or rat serum: Goat anti-Mouse IgG (H&L) conjugated to HRP, affinity purified, min x w/ bovine, horse, human, pig or rabbit serum proteins
Cat. No.:	2UH-100
Amount:	100 mL
Concentration:	Ready to Use
Conjugate:	HRP
Applications:	IHC-P – optimized for use with autostainer or manual staining
Form:	Liquid
Buffer:	TBS, 1% (w/v) BSA, 0.1% Proclin 150, 10% heat-inactivated animal serum
Purification:	Affinity purified using solid phase Goat IgG
Purity:	Affinity purified antibody is > 95% based on SDS-PAGE
Host:	Goat
Clonality:	Polyclonal
Immunogen:	Purified Goat IgG, whole molecule
Storage:	2-8 °C
Shelf-Life:	12 months from received date
Specificity:	Based on IEP, the anti-Rabbit immunopolymer antibody reacts with heavy (γ) chains on rabbit IgG, light chains on all rabbit immunoglobulins: the anti-Mouse immunopolymer antibody reacts with heavy (γ) chains on mouse IgG, light chains on all mouse immunoglobulins
Cross-Reactivity:	Based on IEP, the anti-Rabbit immunopolymer has no cross-reactivity to non-immunoglobulin rabbit serum proteins, serum proteins from human, mouse or rat: the anti-Mouse immunopolymer has no cross-reactivity to non-immunoglobulin mouse serum proteins, serum proteins from bovine, horse, human, pig/swine or rabbit
Country of Origin:	Goat serum was obtained from healthy animals of US origin and under the care of a registered veterinarian.
Disclaimer:	<i>For In vitro Laboratory Use Only. Not for diagnostic or therapeutic use. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Cell IDx, inc. Product may not be resold or modified for resale without prior written approval of Cell IDx, inc. Liability of Cell IDx is limited to the purchase price of the product</i>

Tonsil stained with rabbit anti-CD3 (1ug/ml) and detected with UltraPolymer Universal goat anti-mouse/rabbit-HRP

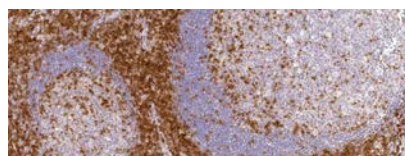

Tonsil stained with mouse anti-Ki-67 (1ug/ml) and detected with UltraPolymer Universal goat anti-mouse/rabbit-HRP

