

---

## Product Specifications


### Cell Palette DAB HRP Chromogen Kit

Description: DAB chromogen for use with HRP enzyme labeled IHC detection systems  
Cat. No.: CPHD-050  
Amount: 50 mL  
Application: IHC

#### Product Specifications

Components: Cell Palette DAB HRP chromogen (50X) and substrate buffer (1X).  
Storage: Chromogen – store at 4°C  
Buffer – for maximum shelf-life store at 4°C. Allow buffer to warm to room temperature prior to use.  
Shelf-Life: 12 months at 4°C  
Instructions: For each 1.0 mL chromogen working solution add 20 µL Chromogen to 980 µL Buffer. Once prepared, the DAB chromogen working solution is stable for 6 hours. Incubate slides for five minutes or until desired intensity has been achieved.

Image: **CD3 (Brown)** staining of FFPE Human Tonsil Tissue

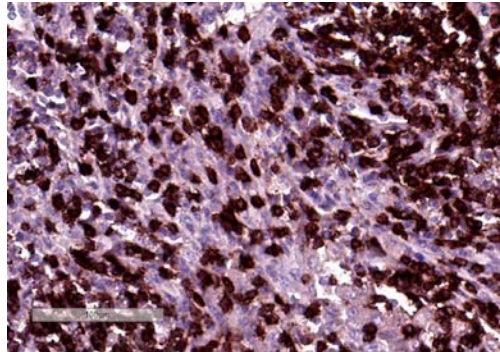

Disclaimer: *For In vitro Research Use Only. Not for diagnostic or therapeutic use. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Cell IDx, Inc. Product may not be resold or modified for resale without prior written approval of Cell IDx, Inc.*

Version 1.0 Date 6.23.21